JAPA tecno

START

CUSTOMIZE

EXIT

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The interface displays a main form for patient data, a search table, a 'New Patient' dialog box, and a dog diagram with measurement points.

Main Form Fields:

- DOG NAME: pp
- OWNER'S NAME: Bob
- AGE: 3 GENDER: FEMALE
- ADDRESS:
- PHONE:
- EMAIL:

Buttons:

- NEW CARD (highlighted with a red box)
- SELECT CARD
- DELETE CARD
- PRESENT ANALYSIS
- EXIT
- RESEARCH
- VIEW RESULTS
- COMPARATIVE ANALYSIS

Search Table:

DATE TIME	PICTURES	NAME
2022/6/18 13:38:00		BODY
2022/6/18 13:38:07		INTESTINE
2022/6/18 13:38:14		LIVER; back view
2022/6/18 13:38:22		BLOOD CELLS
2022/6/18 13:38:29		CRYSTALLINE LENS OF THE EYE
2022/6/18 13:38:34		TRACHEA AND BRONCHI
2022/6/18 13:38:40		LONGITUDINAL SECTION OF RIGHT KIDNEY

New Patient Dialog Box:


DOG NAME: pp
OWNER'S NAME: Bob
GENDER: ☐ MALE ☒ FEMALE
AGE: 3
ADDRESS:
PHONE:
E-MAIL:
OK CANCEL

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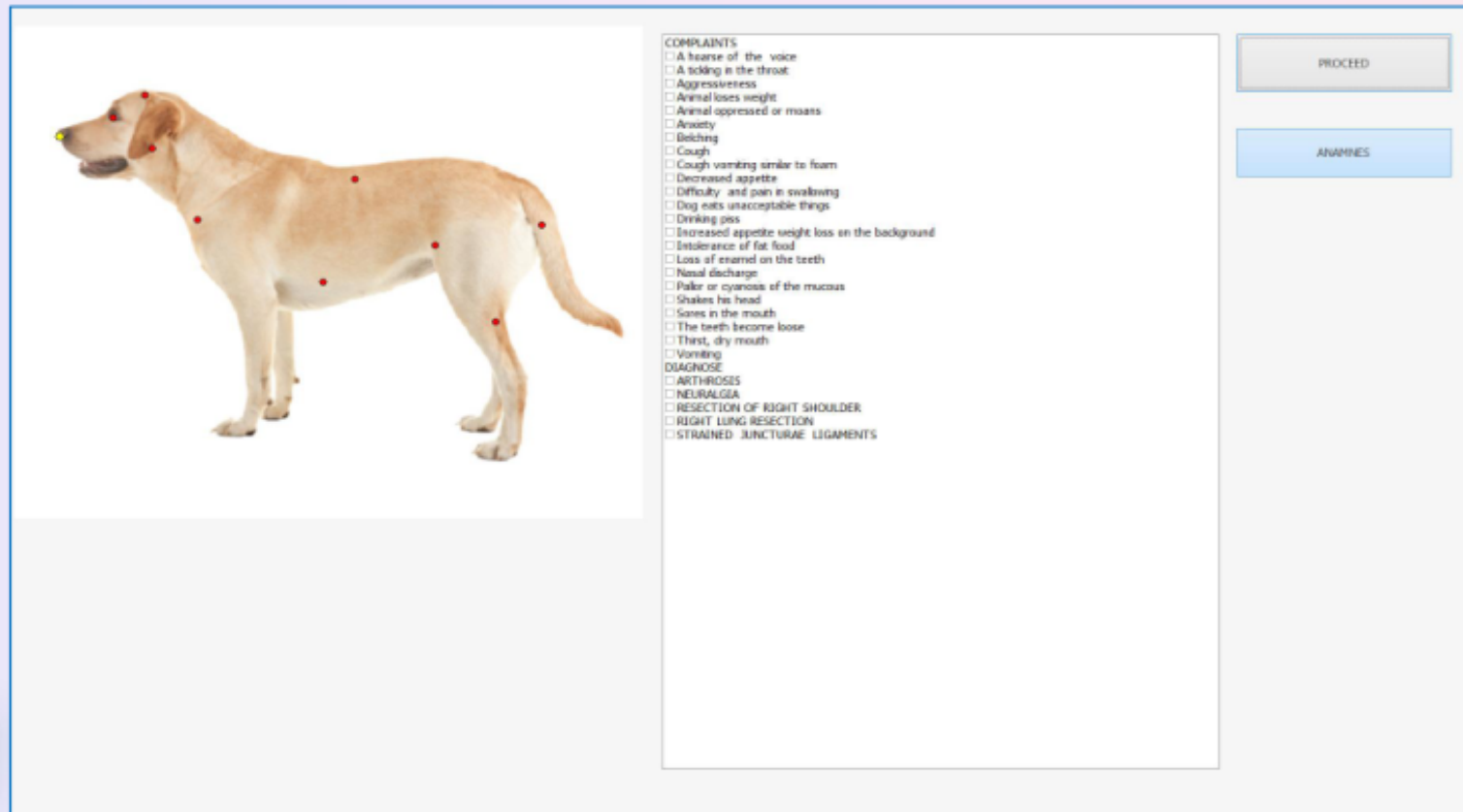COMPLAINTS

- ☐ A hoarse of the voice
- ☐ A tickling in the throat
- ☐ Aggressiveness
- ☐ Animal loses weight
- ☐ Animal oppressed or moans
- ☐ Anxiety
- ☐ Belching
- ☐ Cough
- ☐ Cough vomiting similar to foam
- ☐ Decreased appetite
- ☐ Difficulty and pain in swallowing
- ☐ Dog eats unacceptable things
- ☐ Drinking piss
- ☐ Increased appetite weight loss on the background
- ☐ Intolerance of fat food
- ☐ Loss of enamel on the teeth
- ☐ Nasal discharge
- ☐ Paler or cyanosis of the mucous
- ☐ Shakes his head
- ☐ Sores in the mouth
- ☐ The teeth become loose
- ☐ Thirst, dry mouth
- ☐ Vomiting

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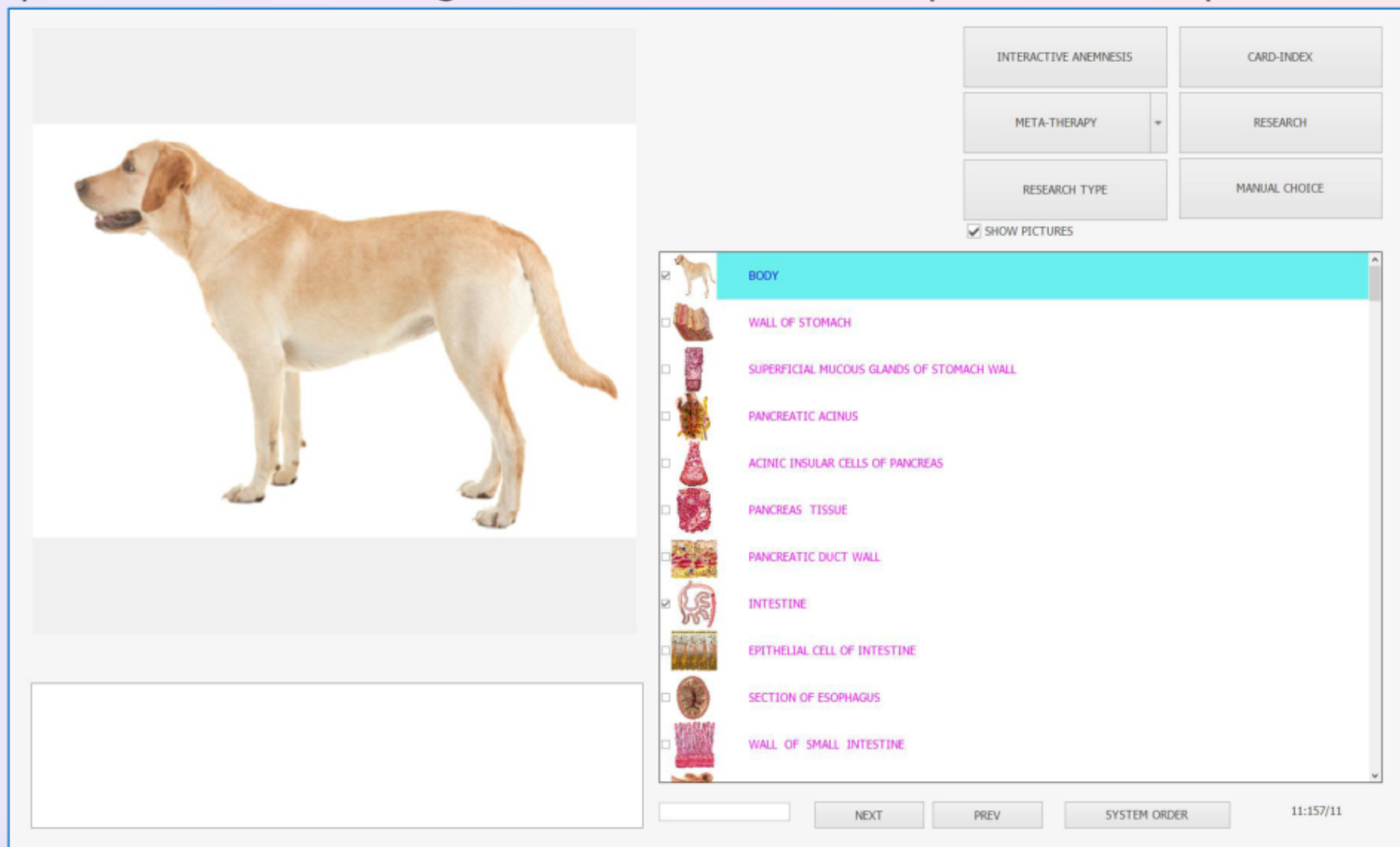
The screenshot displays the NLS diagnostic and therapy software interface. On the left, there is a large image of a dog. Below the image is a text input field. On the right, there is a list of organs and tissues, each with a checkbox and a small icon. The list includes:

- ☒ BODY
- ☐ WALL OF STOMACH
- ☐ SUPERFICIAL MUCOUS GLANDS OF STOMACH WALL
- ☐ PANCREATIC ACINUS
- ☐ ACINIC INSULAR CELLS OF PANCREAS
- ☐ PANCREAS TISSUE
- ☐ PANCREATIC DUCT WALL
- ☒ INTESTINE
- ☐ EPITHELIAL CELL OF INTESTINE
- ☐ SECTION OF ESOPHAGUS
- ☐ WALL OF SMALL INTESTINE

At the top right, there are several buttons: INTERACTIVE ANEMNESIS, CARD-INDEX, META-THERAPY, RESEARCH, RESEARCH TYPE, and MANUAL CHOICE. Below these buttons is a checkbox labeled 'SHOW PICTURES' which is checked. At the bottom, there are buttons for NEXT, PREV, and SYSTEM ORDER, along with a timestamp 11:157/11.

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