JAPA tecno

JAPA tecno

MANUAL DE TREINAMENTO

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* DATE TIME	PICTURES	NAME	ACTIVE	ENERGY	OPTIMUM DISPERSION
2022/10/20 16:56:36		CEREBRUM, top view			1.481
2022/10/20 16:54:20		CHAKRAS			1.988
2022/10/20 16:58:21		CHROMOSOMES (FEMALE)			0.538
2022/10/20 16:54:56		CORONAL CROSS-SECTION OF LARYNX AND TRACHEA			1.605
2022/10/20 16:53:42		CORONAL THORACOTOMY AT THE LEVEL OF ASCENDING PART OF AORTA, FRONT VIEW			2.506
2022/10/21 12:10:28		CORONAL THORACOTOMY AT THE LEVEL OF ASCENDING PART OF AORTA, FRONT VIEW			2.519
2022/10/20 16:53:58		CROSS SECTION OF ABDOMEN AT THE LEVEL OF 1ST LUMBAR VERTEBRA			2.773
2022/10/20 16:58:09		EYE, CRYSTALLINE LENSE			0.791
2022/10/20 16:53:18		HEAD LONGITUDINAL CROSS-SECTION, LEFT			2.431
2022/10/20 16:58:16		HISTAMINE			0.34

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CROSS-SECTION, RIGHT

CROSS-SECTION, LEFT

SECTION OF HEAD AT THE LEVEL OF AQUEDUCT OF CEREBRUM

☐ HORIZONTAL CROSS-SECTION OF HEAD AT THE LEVEL OF THE FOURTH VENTRICLE

☐ HEAD FRONTAL CROSS-SECTION

☐ CROSS - SECTION OF NECK

☒ CORONAL THORACOTOMY AT THE LEVEL OF ASCENDING PART OF AORTA, FRONT VIEW

☐ CORONAL SECTION AT A LEVEL HOLLOW VEIN, FRONT VIEW

☐ HORIZONTAL CROSS-SECTION OF TRUNK AT THE LEVEL OF SHOULDER JOINTS

☐ HORIZONTAL CROSS-SECTION OF CHEST AT THE LEVEL OF 4TH CERVICAL VERTEBRA

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The screenshot shows a software interface with the following components:


- Buttons:** ETALONS, CLEAR, A +, A -
- Category List (Left):**
 - ☐ A # ORGANOPREPARAT
 - ☐ B # BIOCHEMICAL HOMEOSTASIS
 - ☐ C # PATHOLOGY
 - ☐ D # ALLERGENS
 - ☒ E # MICROORGANISMS AND HELMINTHS
 - ☐ F # ALLOPATHY
 - ☐ G # HOMEOPATHY
 - ☐ H # PHYTOTHERAPY
 - ☐ I # NUTRICEVITICS AND PARAPHARMACEUTICALS
 - ☐ J # FOOD
 - ☐ K # LITHOTHERAPY
 - ☐ L # TRACE ELEMENTS
 - ☐ M # THERAPEUTIC COSMETOLOGY
 - ☐ N # BACH FLOWERS
 - ☐ O # PHYSIOTHERAPY
 - ☐ P # EMOTION MATRIX
 - ☐ Q # VITAMIN TABLE
 - ☐ R # HEAVY METALS
 - ☐ S7 SAHASRARA meridians : VC + VG
 - ☐ S6 AJNA meridians : MC + TR
 - ☐ S5 VISHUDHA meridians : P + GI + Sk + Sd
 - ☐ S4 ANAHATA meridians : C + IG
 - ☐ S3 MANIPURA meridians : E + VR + Ed + Ad
- Central Image:** Two purple, rod-shaped bacteria with flagella. A red box highlights this image, and a red arrow points from it to the 'HELICOBACTER PYLORI' entry in the table below.
- Line Graph (Right):** A graph with a black background and green grid lines. The x-axis is labeled with values 1, 2.6, 3.4, 4.2, 4.9, 5.8, 6.6, 7, 8.5. The y-axis has two scales: the left scale ranges from 1 to 6, and the right scale ranges from 8.5 to 240.0. Three lines (red, blue, and yellow) are plotted on the graph.
- Table (Bottom):** A table titled 'ETAON LIST ACCORDING TO DECREASING SPECTRAL SIMILARITY'. It has three columns: a checkbox, a numerical value, and a text description. The entries are:

☐	0 0.000	STOMACH BACK
☐	0 1.439	OPTIMUM DISTRIBUTION
☒	0	Virtual model
☐	0 0.247	HELICOBACTER PYLORI
☐	0 0.412	LACTOBACILLUS ACIDOPHILUS
☐	0 0.728	MUCOR RACEMOSUS

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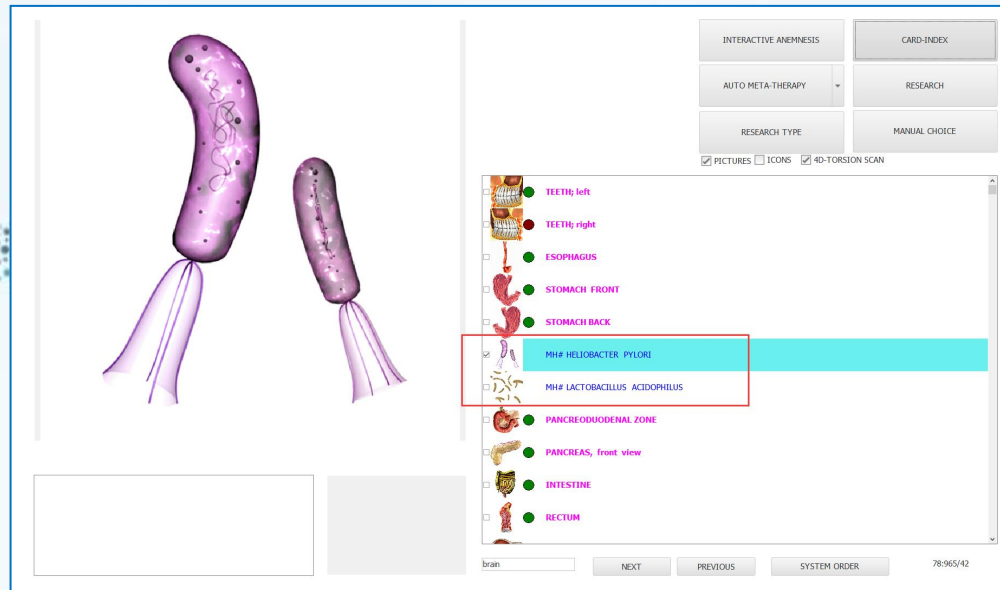
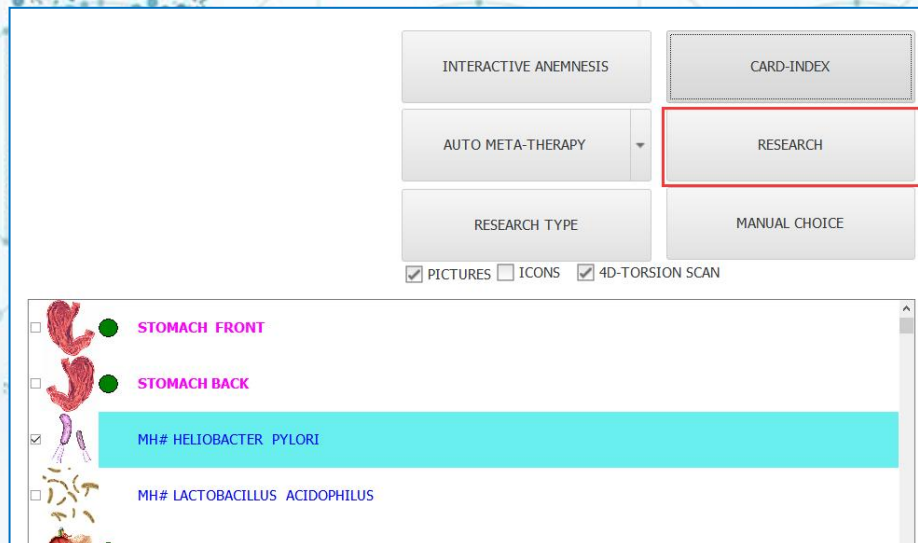
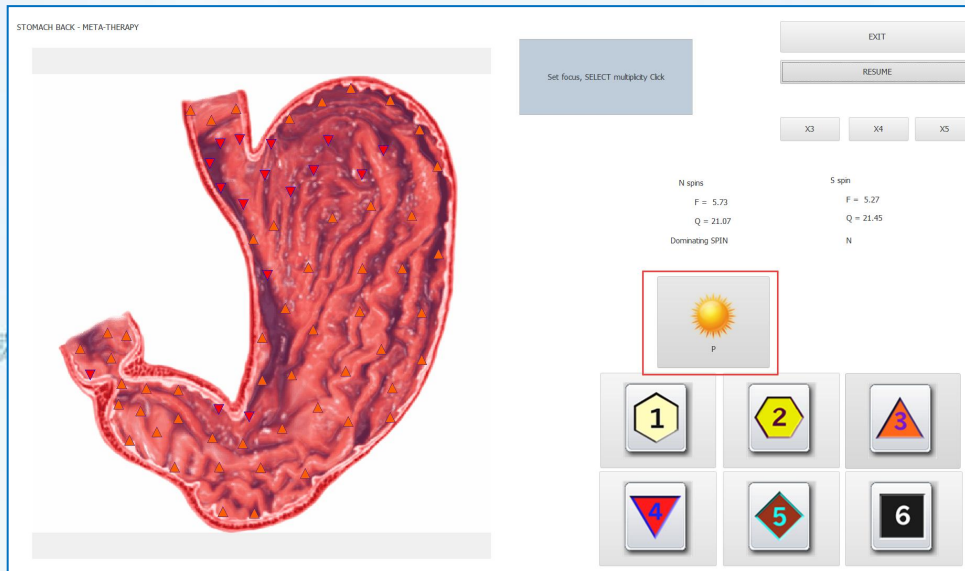
COMPARATIVE ANALYSIS


Enter text to search...FindClear

DATE TIME	PICTURES	NAME	ACTIVE	ENERGY	OPTIMUM DISPERSION
2022/10/24 16:03:02		HEAD LONGITUDINAL CROSS-SECTION, LEFT	S		2,645
2022/10/24 16:03:13		HORIZONTAL CROSS-SECTION OF HEAD AT THE LEVEL OF AQUEDUCT OF CEREBRUM	S		2,097
2022/10/24 16:03:22		CORONAL THORACOTOMY AT THE LEVEL OF ASCENDING PART OF AORTA, FRONT VIEW	S		2,443
2022/10/24 16:03:29		HORIZONTAL THORACOTOMY AT THE LEVEL OF THE 6TH THORACAL VERTEBRA	S		2,175
2022/10/24 16:03:35		CROSS SECTION OF ABDOMEN AT THE LEVEL OF 1ST LUMBAR VERTEBRA	S		2,309
2022/10/24 16:03:40		HORIZONTAL CROSS-SECTION OF TRUNK AT THE LEVEL OF UMBILICUS	S		2,113
2022/10/24 16:03:45		SAGITTAL THORACOTOMY	S		3
2022/10/24 16:03:51		ORGANS OF MALE SMALL PELVIS, right side	S		1,686
2022/10/24 16:03:56		ORGANS OF MALE SMALL PELVIS, left side	S		1,842
2022/10/24 16:04:01		CHAKRAS	S		1,905

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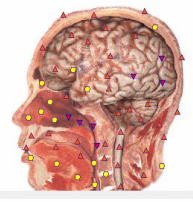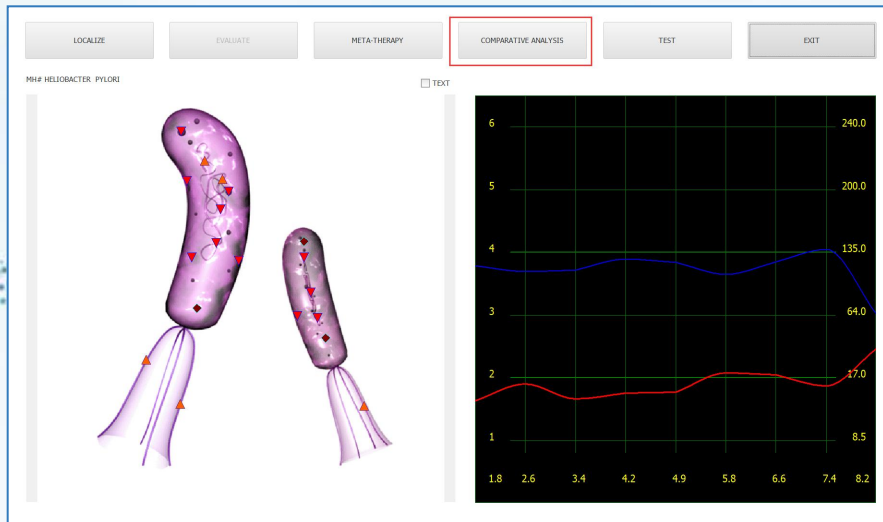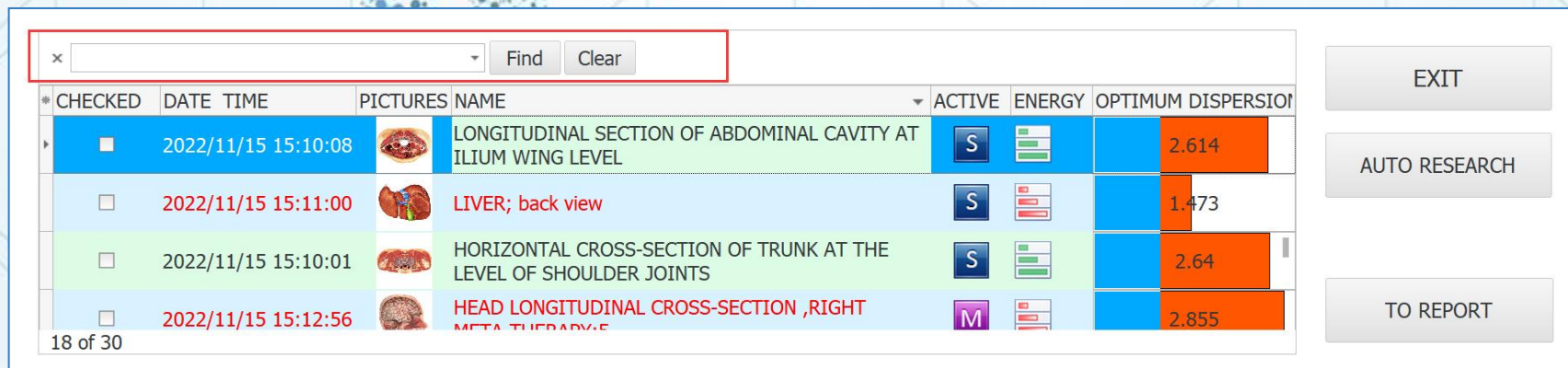
The screenshot displays a medical analysis software interface. At the top, there is a search bar with a magnifying glass icon, a dropdown arrow, and buttons for 'Find' and 'Clear'. Below this is a table with the following columns: *CHECKED, DATE TIME, PICTURES, NAME, ACTIVE, ENERGY, and OPTIMUM DISPERSION. The table contains four rows of data, each with a corresponding medical image in the 'PICTURES' column. The first row is highlighted in blue and has a checked checkbox. The second row is highlighted in light blue and has an unchecked checkbox. The third row is highlighted in light green and has an unchecked checkbox. The fourth row is highlighted in light blue and has an unchecked checkbox. To the right of the table, there are three buttons: 'EXIT', 'AUTO RESEARCH', and 'TO REPORT'. The bottom left corner of the table shows '18 of 30'.

*CHECKED	DATE TIME	PICTURES	NAME	ACTIVE	ENERGY	OPTIMUM DISPERSION
☒	2022/11/15 15:10:08		LONGITUDINAL SECTION OF ABDOMINAL CAVITY AT ILIUM WING LEVEL	S		2.614
☐	2022/11/15 15:11:00		LIVER; back view	S		1.473
☐	2022/11/15 15:10:01		HORIZONTAL CROSS-SECTION OF TRUNK AT THE LEVEL OF SHOULDER JOINTS	S		2.64
☐	2022/11/15 15:12:56		HEAD LONGITUDINAL CROSS-SECTION ,RIGHT	M		2.855

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The screenshot displays the 'ANÁLISE' (ANALYSIS) page of a software application. On the left, there is a list of organs and systems, each with a checkbox and a color-coded label. The central area features a 3D anatomical model of a human torso with internal organs highlighted in green and blue, alongside a line graph with two data series (blue and red) plotted against a grid. The right side of the interface contains a vertical list of analysis options, including 'EXIT', 'V-TEST EXTERNAL', 'V-TEST INTERNAL', 'REPRINTER', 'AURA', 'NLS - ANALYSIS', 'ENTROPY ANALYSIS', 'AUTO MODEL', 'ABSOLUTE MODEL', 'COMPARATIVE ANALYSIS', 'SAVE TO REPORT', and 'VIEW REPORT'. A dialog box is open in the center, titled 'Input name of tested standart', with a text input field and buttons for 'Edit', 'Clear', 'OK', and 'Cancel'.

ETALONS CLEAR A+ A-

☐ A # ORGANOPREPARAT
☐ B # BIOCHEMICAL HOMEOSTASIS
☐ C # PATHOLOGY
☐ D # ALLERGENS
☐ E # MICROORGANISMS AND HELMINTHS
☐ F # ALLOPATHY
☐ G # HOMEOPATHY
☐ H # PHYTOTHERAPY
☐ I # NUTRICEUTICS AND PARAPHARMACEUTICALS
☐ J # FOOD
☐ K # LITHOTHERAPY
☐ L # TRACE ELEMENTS
☐ M # THERAPEUTIC COSMETOLOGY
☐ N # BACH FLOWERS
☐ P # PHYSIOTHERAPY
☐ S7 SAHASRARA meridians : VC + VG
☐ S6 AJNA meridians : MC + TR
☐ S5 VISHUDHA meridians : P + GI + Sk + Sd
☐ S4 ANAHATA meridians : C + IG
☐ S3 MANIPURA meridians : F + VB + Pd + Ad
☐ S2 SVADHISTANA meridians : R + Y + ALL + Ly
☐ S1 MULADHARA meridians : RP + E
☐ VITAMINS INDICATOR
☐ Minerals indicators
☐ EFAs indicators
☐ Toxins indicators
☐ Microbiology indicators
☐ Brain Disorder Table

1.8 SKELETAL SYSTEM
2.6 CONNECTIVE TISSUES, NERVES, CARDIAC VALVES
- LOOSE CONNECTIVE TISSUE, STRIATED MUSCULAR, AND CARDIAC MUSCLE;
3.4 UNSTRATIATED MUSCULAR TISSUE;
4.2 TESSELLATED EPITHELIUM OF THE DIGESTIVE TRACT;
4.9 STRATIFIED SQUAMOUS AND COLUMNAR EPITHELIA.
PARENCHYMATOUS LIVER TISSUE AND TISSUE OF THE BILIARY TRACT;
- KIDNEY TISSUE EPITHELIUM AND REPRODUCTIVE ORGANS;
5.8 LYMPHOID RINGS OF THE PHARYNX, UPPER SECTION OF THE RESPIRATORY TRACT, LYMPHATIC SYSTEM, SPLEEN, OVARIES, AND PROSTATE;
6.6 PERIPHERAL NERVOUS SYSTEM, BRONCHUS EPITHELIUM, ADRENALS, AND THYROID;
7.4 CENTRAL SECTIONS OF SENSORY ANALYSERS EXCEPT THE OPTIC ONES, AND SUB CORTICAL STRUCTURES OF THE BRAIN, PONS CEREBELLI; CEREBELLUM, LIMBIC SYSTEM AND LUNGS PARENCHYMA;
8.2 RETINA, OPTIC NERVE, CEREBRAL CORTEX.

SMART FILTER
SELECT
SHOW ALL
DESCRIPTION
SORT
ADJUSTMENT
ADJUST ALL
AUTO SCAN


Input name of tested standart

Edit Clear OK Cancel

EXIT
V-TEST EXTERNAL
V-TEST INTERNAL
REPRINTER
AURA
NLS - ANALYSIS
ENTROPY ANALYSIS
AUTO MODEL
ABSOLUTE MODEL
COMPARATIVE ANALYSIS
SAVE TO REPORT
VIEW REPORT

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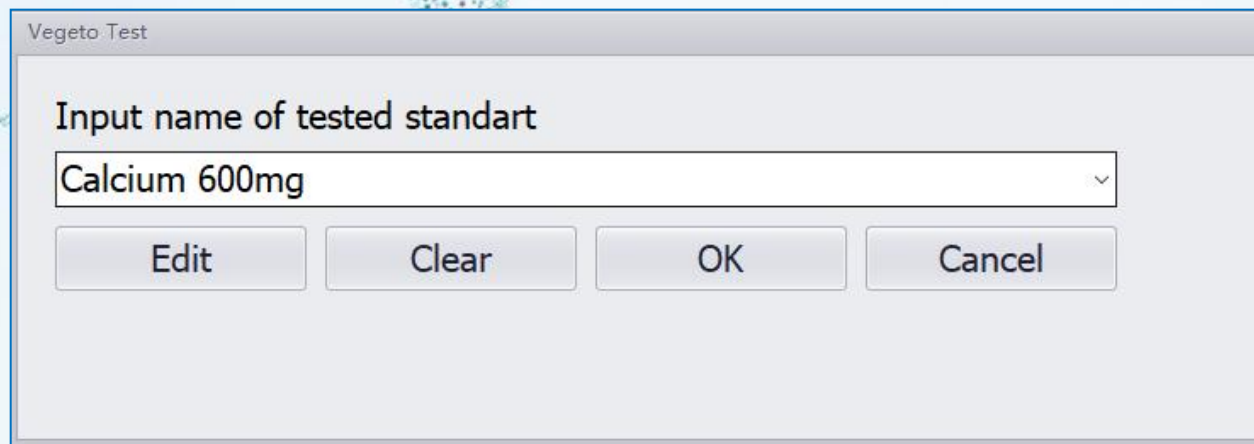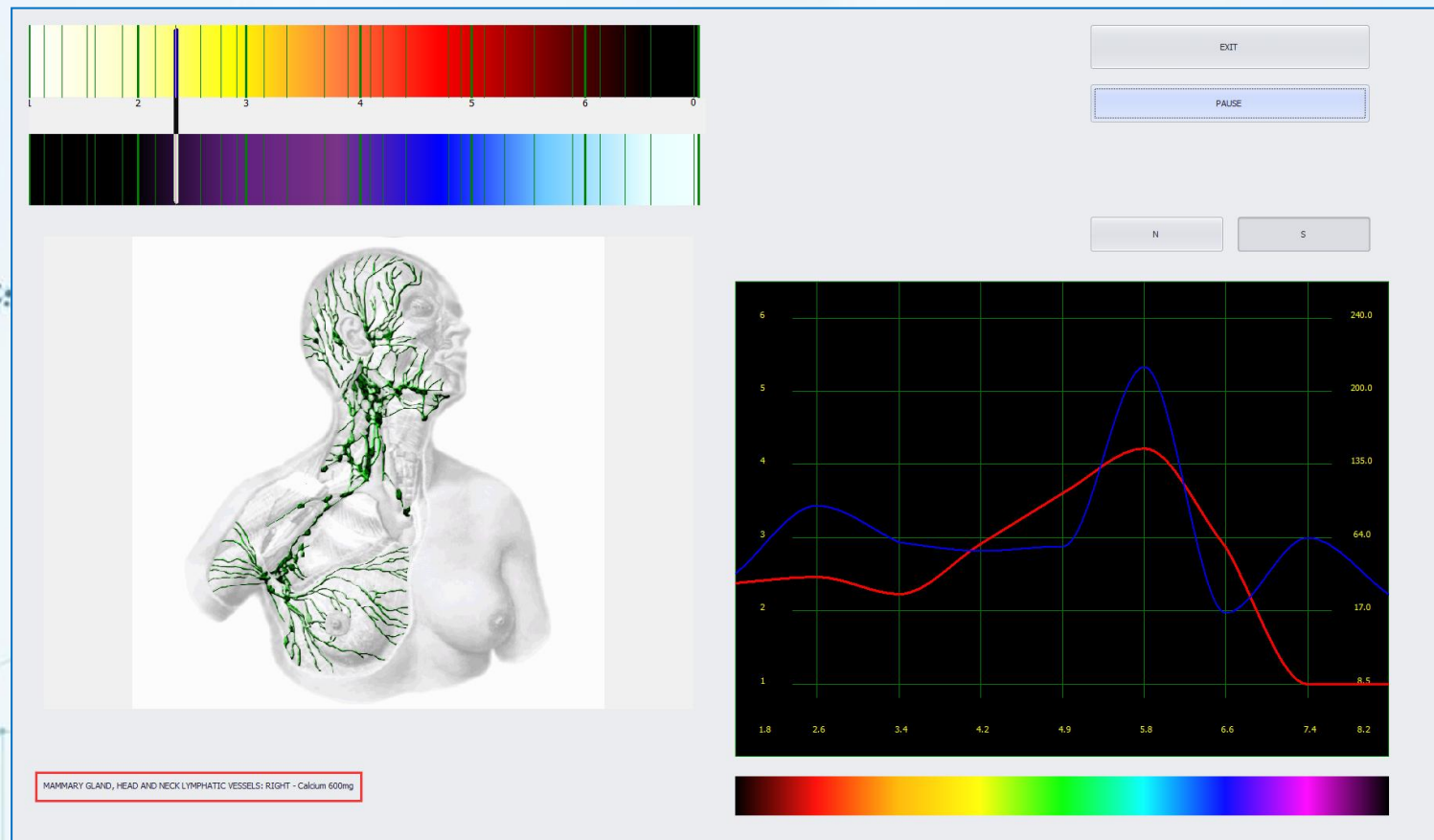
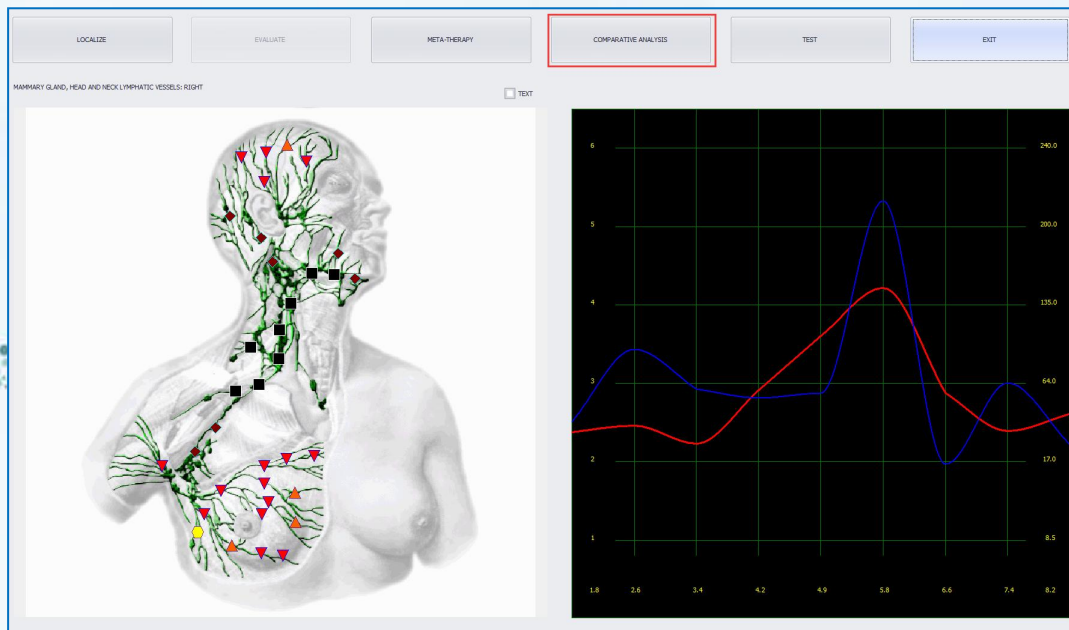
☐ SMART FILTER
☐ SELECT

SHOW ALL DESCRIPTION SORT ADJUSTMENT ADJUST ALL AUTO SCAN

ETALON LIST ACCORDING TO DECREASING SPECTRAL SIMILARITY

0	0.000	MAMMARY GLAND, HEAD AND NECK LYMPHATIC VESSELS: RIGHT
0	1.011	OPTIMAL DISTRIBUTION
X 0		Virtual model

1.8 SKELTAL SYSTEM
2.8 CONNECTIVE TISSUES, NERVES, CARDIAC VALVES
3.8 LOOSE CONNECTIVE TISSUE, STRIATED MUSCULAR, AND CARDIAC MUSCLES
4.8 STRATIFIED MUSCULAR TISSUE
5.8 STRATIFIED EPITHELIUM OF THE DIGESTIVE TRACT
6.8 STRATIFIED SQUAMOUS AND COLUMNAR EPITHELIA, ENDOCRINE GLANDS, LIVER, TISSUE AND TISSUE OF THE BILINARY TRACT
7.8 TRANSITORY EPITHELIUM AND REPRODUCTIVE ORGANS
8.8 LUTHERED RING OF THE PHARYNX, UPPER SECTION OF THE RESPIRATORY TRACT, LYMPHATIC SYSTEM, SPLEEN, CHANGEL, AND PROSTATE
9.8 PERIPHERAL NERVOUS SYSTEM, BRONCHUS EPITHELIUM, ADRENALS, AND THYROID
0.8 CENTRAL SECTIONS OF SENSORY ANALYSES EXCEPT THE OPTIC ONES, AND SUB CORTICAL STRUCTURES OF THE BRAIN, PONS, CEREBELLUM, CEREBELLUM, LIMBIC SYSTEM AND LUNGS, PARSCHYAL
1.8 RETINAL, OPTIC NERVE, CEREBRAL CORTEX

EXIT
V-TEST EXTERNAL
V-TEST INTERNAL
REPORTER
AURA
NLS - ANALYSIS
ENTHROPY ANALYSIS
AUTO MODEL
ABSOLUTE MODEL
COMPARATIVE ANALYSIS
SAVE TO REPORT
VIEW REPORT

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0	0.000	MAMMARY GLAND, HEAD AND NECK LYMPHATIC VESSELS: LEFT
0	0.886	OPTIMUM DISTRIBUTION
x	0	Virtual model
0	0.422	MONONUCLEOSIS
0	0.678	LYMPHANGITIS
0	0.774	SYSTEMIC LUPUS ERYTHEMATOSUS
0	0.784	SCLERODERMA
0	0.852	LYMPHADENITIS
0	0.882	LIPOMA OF THE MAMMARY GLAND

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The screenshot displays the AURA software interface. On the left, there is a list of organs and chakras with checkboxes. The central part shows a human figure with chakras and a spectral graph. The right side contains a list of organs and chakras, with the 'AURA' button highlighted in red. Below the graph, there is a table with columns for 'CHAKRAS' and 'OPTIMUM DISTRIBUTION'.

CHAKRAS	OPTIMUM DISTRIBUTION
0	0.000
0	1.923
X	0

Below the table, there is a section titled 'ETADON LIST ACCORDING TO DECREASING SPECTRAL SIMILARITY'.

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
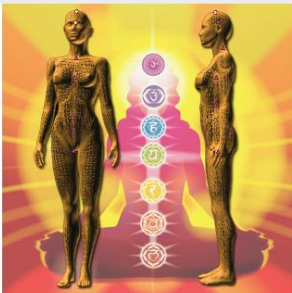
ADJUST ALL

AUTO SCAN

ETAON LIST ACCORDING TO DECREASING SPECTRAL SIMILARITY

0	0.000	CHAKRAS
0	1.923	OPTIMUM DISTRIBUTION
X 0		Virtual model

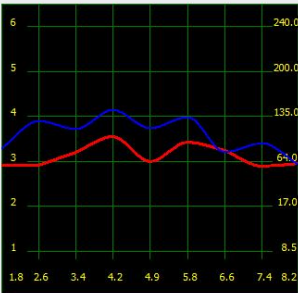1.8 SKELETAL SYSTEM
2.6 CONNECTIVE TISSUES, NERVES, CARDIAC VALVES
- LOOSE CONNECTIVE TISSUE, STRIATED MUSCULAR, AND CARDIAC MUSCLE;
3.4 UNSTRIATED MUSCULAR TISSUE.
4.2 TESSELLATED EPITHELIUM OF THE DIGESTIVE TRACT;
4.9 STRATIFIED SQUAMOUS AND COLUMNAR EPITHELIA, PARENCHYMATOUS LIVER TISSUE AND TISSUE OF THE BILIARY TRACT;
- KIDNEY TISSUE EPITHELIUM AND REPRODUCTIVE ORGANS;
5.8 LYMPHOID RING OF THE PHARYNX, UPPER SECTION OF THE RESPIRATORY TRACT, LYMPHATIC SYSTEM, SPLEEN, OVARIES, AND PROSTATE;
6.6 PERIPHERAL NERVOUS SYSTEM, BRONCHUS EPITHELIUM, ADRENALS, AND THYROID;
7.4 CENTRAL SECTIONS OF SENSORY ANALYSERS EXCEPT THE OPTIC ONES, AND SUB CORTICAL STRUCTURES OF THE BRAIN, PONS CEREBELLI, CEREBELLUM, LIMBIC SYSTEM AND LUNGS PARENCHYMA;
8.2 RETINA, OPTIC NERVE, CEREBRAL CORTEX.

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The screenshot displays the JAPA tecno software interface for the Fitoterapia (Phytotherapy) section. On the left, a list of therapy options is shown, with 'H # PHYTOTHERAPY' selected and highlighted. The central area features a photograph of a plant, identified as Oxalis acetosella. To the right of the plant image is a spectral graph showing two curves (red and blue) plotted against frequency (1.8 to 8.2) and amplitude (0 to 6). Below the graph, a table titled 'ETAON LIST ACCORDING TO DECREASING SPECTRAL SIMILARITY' lists various plant species and their corresponding coefficients. The table is as follows:

0	0.000	NERVES OF BACK
0	1.325	OPTIMUM DISTRIBUTION
X 0		Virtual model
0	0.853	OXALIS ACETOSELLA*
0	0.868	RHODIOLA ROSEA*
0	0.891	ASPARAGUS OFFICINALIS*
0	0.911	RIBUS NIGRUM*
0	0.920	VERBASCUM THAPSUS*
0	0.937	CORYLUS AVELLANA*

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